

For Models:

**DCRV75**

Vaccine refrigerators

**1** Unpack the refrigerator, solar panels, battery box, and wiring materials. (make sure no damage incurred in shipping).



**2** Safely install solar panels on roof using supplied framing brackets. Make sure the panels are not shaded.

**3** Place refrigerator in a cool location, out of the sunlight, keeping adequate ventilation space between the refrigerator and walls.

**4** Attach battery box on the wall near the refrigerator in secure manner.

**5** Using supplied wires, attach solar array to the battery box and the battery box to the refrigerator as shown below. (The solar panels are connected in parallel, not series).

**6** Make sure switch is turned on. Wait 2 days before loading refrigerator with vaccines to give time for the refrigerator to get cold and freeze the interior packs.

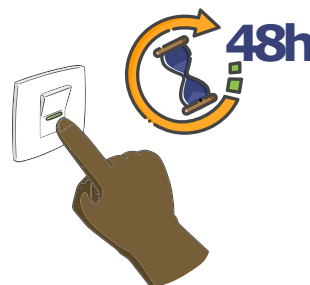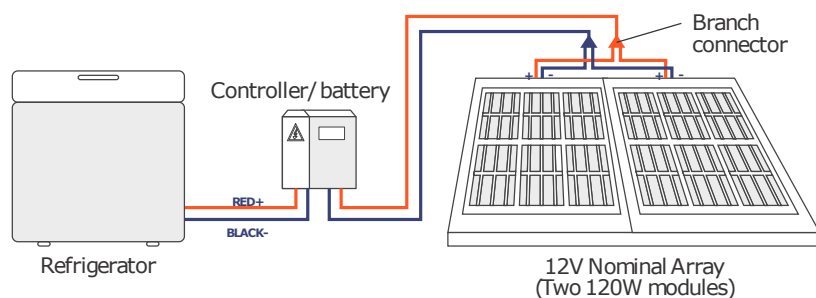

Exclusively distributed by :  
**ElitesGroup Inc.**

651 N. Broad St., Suite 206, Middletown, DE 19709  
Cell : +1 (510) 306 2359  
Email: [contact@elitesgroup.co](mailto:contact@elitesgroup.co)  
Web : [www.elitesgroup.co](http://www.elitesgroup.co)

**elitesGroup**